

# Norovirus VP1 antibody [HL2177]

## Cat. No. GTX638181

Host	Rabbit
Clonality	Monoclonal
Isotype	IgG
Application	WB, ICC/IF
Reactivity	Norovirus

Package  
100 µl, 25 µl

### APPLICATION

#### Application Note

\*Optimal dilutions/concentrations should be determined by the researcher.

Suggested dilution	Recommended dilution
WB	1:500-1:3000
ICC/IF	Assay dependent

Not tested in other applications.

**Product Note** This antibody was raised against Norovirus GII.17 VP1 and does not cross-react with Norovirus GII.4 VP1.

### PROPERTIES

Form	Liquid
Buffer	PBS
Preservative	No preservative
Storage	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
Concentration	2 mg/ml (Please refer to the vial label for the specific concentration.)
Immunogen	Synthetic peptide of Norovirus VP1. (Norovirus Hu/GII/CN/2014/GII.17/1401169)
Purification	Affinity purified by Protein A.
Conjugation	Unconjugated

**Note** For laboratory research use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.

Purchasers shall not, and agree not to enable third parties to, analyze, copy, reverse engineer or otherwise attempt to determine the structure or sequence of the product.

Scan me for more info!



For full product information, images and publications, please visit our [website](https://www.genetex.com).



[www.unibiotech.in](https://www.unibiotech.in)



An ISO 9001:2015 Certified Company


Order now and enjoy our services anywhere in India!

Tel: +91 9910 73 8844 | Email: [info@unibiotech.in](mailto:info@unibiotech.in)

Reg. Address : A-33, 2nd Floor, Ring Road, Rajouri Garden, New Delhi-110027 (INDIA)


Date 2024 / 12 / 18 Page 1 of 2

DATA IMAGES

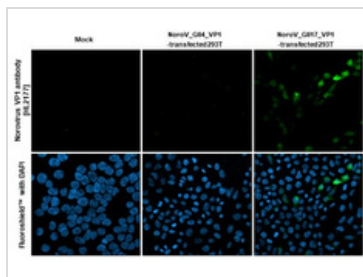

GTX638181 ICC/IF Image

Norovirus VP1 antibody [HL2177] detects Norovirus VP1 protein by immunofluorescent analysis. Sample: Mock and transfected 293T cells were fixed in 4% paraformaldehyde at RT for 15 min. Green: Norovirus VP1 stained by Norovirus VP1 antibody [HL2177] (GTX638181) diluted at 1:500. Blue: Fluoroshield with DAPI (GTX30920).

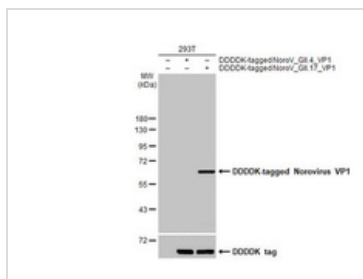

GTX638181 WB Image

Non-transfected (-) and transfected (+) 293T whole cell extracts (30 µg) were separated by 10% SDS-PAGE, and the membrane was blotted with Norovirus VP1 antibody [HL2177] (GTX638181) diluted at 1:2000. The HRP-conjugated anti-rabbit IgG antibody (GTX213110-01) was used to detect the primary antibody.

Scan me for more info!



For full product information, images and publications, please visit our [website](http://www.unibiotech.in).



[www.unibiotech.in](http://www.unibiotech.in)



An ISO 9001:2015 Certified Company

Order now and enjoy our services anywhere in India!

Tel: +91 9910 73 8844 | Email: [info@unibiotech.in](mailto:info@unibiotech.in)

Reg. Address : A-33, 2nd Floor, Ring Road, Rajouri Garden, New Delhi-110027 (INDIA)

Date 2024 / 12 / 18 Page 2 of 2