

Somatostatin antibody [HMV303] HistoMAX™

Cat. No. GTX639926

Host	Rabbit
Clonality	Monoclonal
Isotype	IgG
Application	IHC-P
Reactivity	Human

Package

500 µl, 100 µl

PRODUCT

This antibody was validated on 76 different Normal Tissues by IHC-P.

Summary

[Go to Normal Tissue Gallery](#)
[Go to Cancer Tissue Gallery](#)

APPLICATION

Application Note

*Optimal dilutions/concentrations should be determined by the researcher.

Suggested dilution	Recommended dilution
IHC-P	1:100-1:200

Note : Heat-induced antigen retrieval for 5 minutes in an autoclave at 121°C in pH 7.8 Target Retrieval Solution buffer.

Not tested in other applications.

PROPERTIES

Form	Liquid
Buffer	PBS, 2% BSA
Preservative	0.025% ProClin 300
Storage	Store as concentrated solution. Centrifuge briefly prior to opening vial. For short-term storage (1-2 weeks), store at 4°C. For long-term storage, aliquot and store at -20°C or below. Avoid multiple freeze-thaw cycles.
Concentration	28 µg/ml (Please refer to the vial label for the specific concentration.)
Immunogen	Carrier-protein conjugated synthetic peptide encompassing a sequence within the center region of mouse Somatostatin. The exact sequence is proprietary.
Purification	Affinity purified by Protein A.
Conjugation	Unconjugated

Scan me for more info!



For full product information, images and publications, please visit our [website](#).



www.unibiotech.in



An ISO 9001:2015 Certified Company

Order now and enjoy our services anywhere in India!

Tel: +91 9910 73 8844 | Email: info@unibiotech.in


Reg. Address : A-33, 2nd Floor, Ring Road, Rajouri Garden, New Delhi-110027 (INDIA)

Date 2024 / 12 / 18 Page 1 of 2

Note


For laboratory research use only. Not for any clinical, therapeutic, or diagnostic use in humans or animals. Not for animal or human consumption.
Purchasers shall not, and agree not to enable third parties to, analyze, copy, reverse engineer or otherwise attempt to determine the structure or sequence of the product.

DATA IMAGES

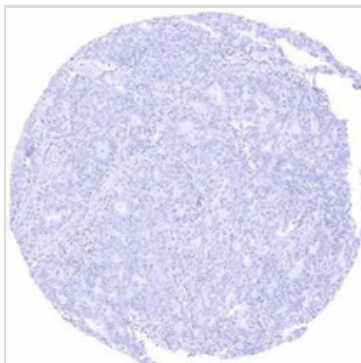


GTX639926 IHC-P Image

IHC-P analysis of human ovarian high-grade serous carcinoma (HGSC) tissue using GTX639926 Somatostatin antibody [HMV303] HistoMAX™.
Somatostatin negative serous high-grade carcinoma.

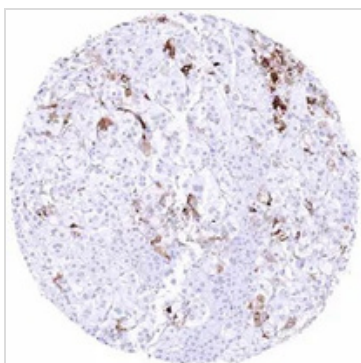

GTX639926 IHC-P Image

IHC-P analysis of human adrenal gland tissue using GTX639926 Somatostatin antibody [HMV303] HistoMAX™.
A rather small subset of cells shows a strong somatostatin positivity.

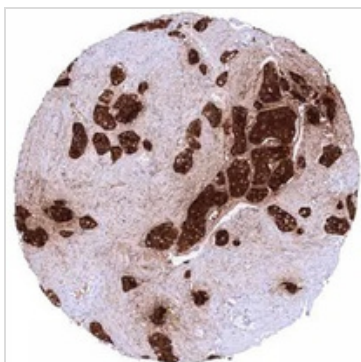

GTX639926 IHC-P Image

IHC-P analysis of human duodenal neuroendocrine tumor (d-NET) tissue using GTX639926 Somatostatin antibody [HMV303] HistoMAX™.
Neuroendocrine tumor with strong somatostatin staining of all tumor cells.



For full product information, images and publications, please visit our [website](https://www.unibiotech.in).

Scan me for more info!



www.unibiotech.in



Order now and enjoy our services anywhere in India!

Tel: +91 9910 73 8844 | Email: info@unibiotech.in

Reg. Address : A-33, 2nd Floor, Ring Road, Rajouri Garden, New Delhi-110027 (INDIA)

Date 2024 / 12 / 18 Page 2 of 2