

Recombinant Human Interleukin-18/IL-18 Protein (His Tag)

Catalog Number: PKSH032625

Note: Centrifuge before opening to ensure complete recovery of vial contents.


Description

Species	Human
Mol_Mass	19.3 kDa
Accession	Q14116
Bio-activity	Not validated for activity

Properties

Purity	> 95 % as determined by reducing SDS-PAGE.
Endotoxin	< 0.01 EU per µg of the protein as determined by the LAL method.
Storage	Generally, lyophilized proteins are stable for up to 12 months when stored at -20 to -80°C. Reconstituted protein solution can be stored at 4-8°C for 2-7 days. Aliquots of reconstituted samples are stable at < -20°C for 3 months.
Shipping	This product is provided as lyophilized powder which is shipped with ice packs.
Formulation	Lyophilized from a 0.2 µm filtered solution of 20mM PB, 150mM NaCl, pH7.4. Normally 5% - 8% trehalose, mannitol and 0.01% Tween 80 are added as protectants before lyophilization. Please refer to the specific buffer information in the printed manual.
Reconstitution	Please refer to the printed manual for detailed information.

Data

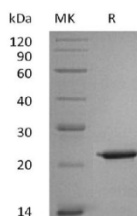

> 95 % as determined by reducing SDS-PAGE.

Background

Interleukin-18 is a secreted protein and it belongs to the IL-1 family. IL-18 is a proinflammatory cytokine and produced by macrophages and other cells. This cytokine can induce the IFN-gamma production of T cells. The combination of this cytokine and IL12 has been shown to inhibit IL-4 dependent IgE and IgG1 production; and enhance IgG2a production of B cells. IL-18 binding protein (IL18BP) can specifically interact with this cytokine; and thus negatively regulate its biological activity. After stimulation with IL-18; natural killer (NK) cells and certain T cells release another important cytokine called interferon- γ (IFN- γ) or type II interferon that plays an important role in activating the macrophages or other cells.

For Research Use Only