

## Purified Anti-Human CD27 Antibody[O323]

**Catalog Number:**E-AB-F11400P

**Note:** Centrifuge before opening to ensure complete recovery of vial contents.

### Description


<b>Reactivity</b>	Human
<b>Immunogen</b>	Recombinant Human CD27 protein
<b>Host</b>	Mouse
<b>Isotype</b>	Mouse IgG1, $\kappa$
<b>Clone</b>	O323
<b>Purification</b>	>98%, Protein A/G purified
<b>Conjugation</b>	Unconjugated
<b>Buffer</b>	PBS, pH 7.2. Contains 0.05% proclin 300.

### Applications

### Recommended Dilution

**FCM** 2  $\mu\text{g/mL}$  ( $1 \times 10^5$ – $5 \times 10^5$  cells)

### Data

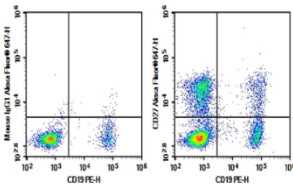

Human peripheral blood lymphocytes were stained with 0.2  $\mu\text{g}$  Purified Anti-Human CD27 Antibody[O323] (Right) and 0.2  $\mu\text{g}$  Mouse IgG1,  $\kappa$  Isotype Control (Left), followed by Alexa Fluor® 647-conjugated Goat Anti-Mouse IgG Secondary Antibody, then anti-Human CD19 PE-conjugated Monoclonal Antibody.

### Preparation & Storage

<b>Storage</b>	Store at 4°C valid for 12 months or -20°C valid for long term storage, avoid freeze / thaw
<b>Shipping</b>	Ice bag

### Background

The protein encoded by this gene is a member of the TNF-receptor superfamily. This receptor is required for generation and long-term maintenance of T cell immunity. It binds to ligand CD70, and plays a key role in regulating B-cell activation and immunoglobulin synthesis. This receptor transduces signals that lead to the activation of NF-kappaB and MAPK8/JNK. Adaptor proteins TRAF2 and TRAF5 have been shown to mediate the signaling process of this receptor. CD27-binding protein (SIVA), a proapoptotic protein, can bind to this receptor and is thought to play an important role in the apoptosis induced by this receptor.

### For Research Use Only