

Recombinant Human TCbIR/8D6A Protein (His Tag)

Catalog Number: PKSH033133

Note: Centrifuge before opening to ensure complete recovery of vial contents.


Description

Species	Human
Mol_Mass	21.1 kDa
Accession	Q9NPF0
Bio-activity	Not validated for activity

Properties

Purity	> 95 % as determined by reducing SDS-PAGE.
Endotoxin	< 1.0 EU per µg of the protein as determined by the LAL method.
Storage	Generally, lyophilized proteins are stable for up to 12 months when stored at -20 to -80°C. Reconstituted protein solution can be stored at 4-8°C for 2-7 days. Aliquots of reconstituted samples are stable at < -20°C for 3 months.
Shipping	This product is provided as lyophilized powder which is shipped with ice packs.
Formulation	Lyophilized from a 0.2 µm filtered solution of 10mM Tris-Citrate, 150mM NaCl, pH 7.4. Normally 5% - 8% trehalose, mannitol and 0.01% Tween 80 are added as protectants before lyophilization. Please refer to the specific buffer information in the printed manual.
Reconstitution	Please refer to the printed manual for detailed information.

Data

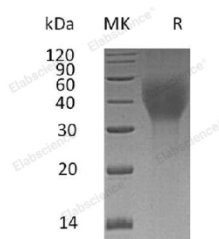

> 95 % as determined by reducing SDS-PAGE.

Background

CD320 antigen is also known as 8D6 antigen, FDC-signaling molecule 8D6, Transcobalamin receptor and 8D6A. It is a single-pass type I membrane protein and containing two LDL-receptor class A domains. CD320 has been recently discovered and reported as a follicular dendritic cell (FDC) protein. CD320 can augment the proliferation of plasma cells precursors generated by IL-10. CD320 also functions as a receptor for the cellular uptake of transcobalamin bound cobalamin. Defects in CD320 are the cause of methylmalonic aciduria type TCbIR (MMATC) which is a metabolic disorder.

For Research Use Only